



Red-Light Therapy Belt

Instruction Manual

FOR INSTRUCTIONS IN DE, FR, IT, ES
VISIT WWW.NEBULA-APPLIANCE.COM

Introduction to you new Red Light Therapy Belt

Before starting your first session, please make sure to read this manual as it explains important operating procedures.

Notes: Each LED inside has three chips at 2pcs 660nm and 1 pcs 850nm
Near Infrared 850nm is invisible to the human eye

This Package Includes:

- ⊙ 1 Red light therapy pad
- ⊙ 1 Controller
- ⊙ 1 Adaptor
- ⊙ 1 Strap
- ⊙ 1 Connecting wire
- ⊙ 1 User manual



Red Light
Therapy Pad



Strap



Adapter



Controller



USB Wire

Specification

Rating Power :	60W
Actual Power:	20W
Product Size:	445*200*4.95MM
LED QTY:	120*0.5W
Net Weight:	0.58KG
Wavelength:	2*660nm; 1*850nm
EMF:	0uT
Timer:	5/10/15/20/25/30 Mins.
Material:	Fabric
Lifespan:	Over 50,000 Hours

The red-light device uses combination red 660nm and infrared 850nm LEDs which are chosen to give the ideal wavelengths found by the scientific research of the last three decades.

660nm Of Red Light

This light can be seen by the human eye and targets outer skin cells

850nm of near infrared-red light

This light cannot be seen by the human eye and targets deep tissue

User Instruction

- 01** Plug into the AC socket Or use the right power bank to power the device. Power bank output required: equal to or exceed 5V 2A.
- 02** **Distance:** 0 inch for best results, place in any area you would like to treat.
- 03** Select Options



Top button: Brightness adjust

Left button: Red/Infrared/Both on

Middle button: 10Hz/40Hz/ Breathe and pulse rate

Right button: 5-30 mins (5mins/unit) timer

Bottom Button: ON/OFF

04 Usage frequency: Use the device daily 1-2 times and for at least 4 weeks, and each time **we** suggest 20-30 minutes on one treatment area. Then move the device to the next area you would like to treat.

05 Unplug the device and store it in a dry area

06 Maintenance: Continue using several times a week

THE RED-LIGHT THERAPY PAD CAN BE USED IN 2 WAYS



2 Different Power Modes

Benefits backed by Science:

Research shows that the wavelengths in the spectrum of 660-900nm will give the optimal biological response.

- ⊙ Muscle and joints pain
- ⊙ Speed healing
- ⊙ Hair Re-growth
- ⊙ Enhanced blood circulation
- ⊙ Aging spots and wrinkle removal
- ⊙ Increase collagen
- ⊙ Weight loss
- ⊙ Anti-inflammatory
- ⊙ Decreases pain

This technology does not have any known side effects but we suggest taking the following precautions:

- PLEASE NOTE RESULTS ARE NOT INSTANT YOU WILL NEED TO USE CONSISTANTLY TO SEE BENEFITS.
- Clear the area being treated by removing all traces of make-up or sweat. For best results wash the area being treated to ensure a clean surface.
- Always unplug the power supply when not in use.
- Do not use device if your skin is light sensitive, if you are currently using medication, or if you're undergoing treatment that may cause sensitivity of skin to light
- Ask your doctor prior to using if you:
 - Are pregnant or suspect you may be pregnant.
 - Are taking drugs known to cause light sensitivity


- Are suffering from a serious medical condition
- Not intended for use with children
- Do not exceed 30-minute treatment
- Stop using immediately and contact us if you suspect that you have had a negative reaction from using this device contact us: - sales@nebula-appliances.com

Portable and flexible

Portable and flexible design you can enjoy it any time, anywhere by connecting a power bank or charging terminal



Application:

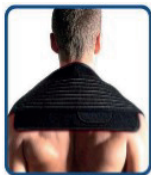


NECK



SHOULDER



KNEE

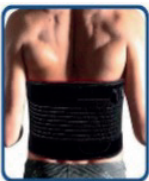

BACK



Relaxing



Working



Meditating