



LED RED LIGHT FACE MASK

NB75



User Manual

www.nebula-appliances.com

LED FACIAL MASK:



The Nebula LED Red Light Face Mask is a wearable and rechargeable beauty instrument for home use that harnesses the power of LED phototherapy. The surface of this product is made of medical-grade silicone, making it convenient, user-friendly, and easily foldable for portability. It can be adjusted to fit different face shapes with its nylon buckle fixation, catering to a variety of users.

Equipped with 240 high-quality LED diodes grouped into 60 blocks, with each block featuring 4 LED beads emitting 4 distinct wavelengths (red light at 630nm, blue light at 460nm, yellow light at 590nm, and near-infrared light at 850nm), this device offers versatility. The controller allows for adjustment of energy levels and light duration, providing users with tailored phototherapy sessions.

MODE EFFICACY:

Red Light: 630nm, anti-aging, promote collagen, tightness, increase elastic
Blue Light: 460nm, kill bacteria, acne bacteria, reduce inflammation and scar, increase skin oxygen
Yellow Light: 590nm, balance skin texture, stimulate red blood cells to reduce redness and fine lines
Near-Infrared Light: 850nm, repair the skin, improve blood circulation, increase skin moisture

Infrared light is independently controlled and can also be used with any kind of light.

FIRST-TIME SETTINGS:

Prior to unboxing the package, make sure to fully charge the Nebula LED Face Mask by connecting it to the charger. The controller's indicator will flash during charging to indicate the progress. Once the charge is complete, the battery display will turn off. Allow 4-5 hours for a complete charge before enjoying the comfortable fit of the mask, which can be adjusted using the velcro buttons. When fitting the Nebula mask, thread the strap through the designated inlet on the opposite side of the mask, secure it around the head, and adjust for a snug fit.



Require assistance? Reach out to Nebula Appliances at sales@nebula-appliances.com for support.


Alternatively, you can get in touch through the contact form on our website.

PRODUCT OPERATION EXPLANATION:



1. Power button: long press for 1-2s to enter the on / off state;
2. E key is the energy adjustment key: short press to adjust and switch 3 gear energy output mode;
3. Light time timing key: press the T + button to increase the light time Press the T-button to reduce the light time 5 gears (10 / 15 / 20 / 20 / 25 / 30min) single use recommended for 15-25 minutes.
4. During charging, the battery display key flashes and the battery display key is off after fully charged.

Power - 48W

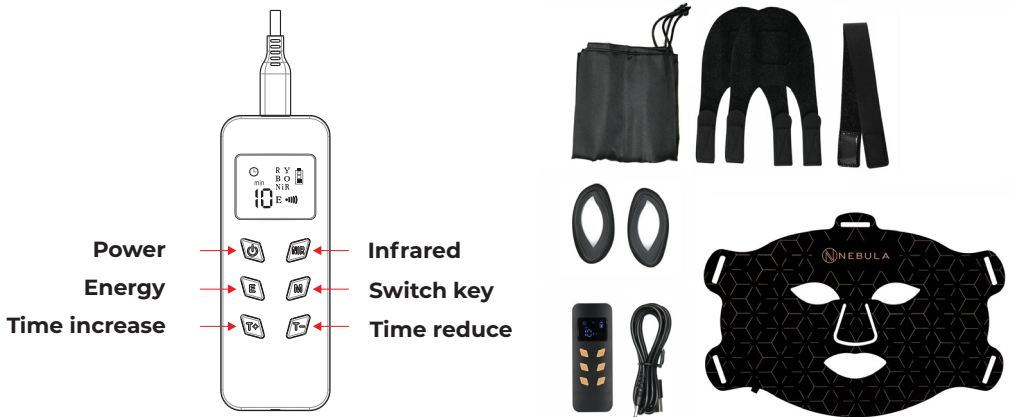

PRODUCT PARAMETER:

Product Name	Nebula LED Red Light Face Mask (NB75)		
LED Lamp Bead	240	External Box Size	43x23x4.5cm
Timer	10/15/20/25/30min	Exchange	0.3A 50/60HZ
Temperature	0-40 C	Net Weight	220g
Adapter Input	100-240V	Adapter Output	5V DC 3000MA
Four Wavelengths Of Light	630nm, 460nm, 590nm, and Near-infrared light 850nm		
Box Contents	Nebula LED Red Light Face Mask, Controller, Charging cable, Instructions Manual, Fixed Strap		

- 1. Please utilise the product in strict accordance with the correct method; do not utilise the equipment outside the product user manual. Doing so could result in machine failure or an accident.**
- 2. Do not operate face mask with damaged cord or plug or after the mask malfunctions, or is dropped or damaged in any manner. Return appliance to the nearest authorised service facility for examination, repair or electrical or mechanical adjustment. Contact Nebula Support for further assistance.**
3. Do not use this product to enhance any situation other than the indication.
4. Pregnant women, those in breastfeeding, and children cannot use this product.
5. If the user is experiencing physical, sensory, or mental illness, or is undergoing dermatological disease improvement, including during the onset period (such as lupus erythematosus, photosensitive eczema, or albinism, etc.), it is advised against using the product. Usage should be strictly under the guidance of a responsible individual, for instance, a skin specialist.
6. If you have any photosensitive disorder, refrain from using this product. Photosensitivity can result in a severe skin reaction.
7. Can not be used on the same day with the photosensitive ingredients.
8. Ensure proper eye protection by correctly wearing the eye mask. It is advisable to keep the eyes closed during use. If you have a genetic eye disease, avoid using a biological mask.
9. If you experience any discomfort while using this product, please cease use immediately. Contact nebula support on **sales@nebula-appliances.com**.
10. This product is not waterproof; do not immerse the equipment in water. If you suspect it has come into contact with water, refrain from using it.
11. Do not give to infants or children; store in childproof areas.
12. Do not use this product if any parts are damaged, including product cabling.
13. This product is intended for daily use only; do not disassemble, modify, or repair without authorisation.
14. We recommend using this product 3-4 times per week, with each use lasting 15-25 minutes.
15. Product Warranty: 1 year standard warranty included with a new device.
Note: For defects or faults, contact our customer support team at **sales@nebula-appliances.com / sales@junglehutonline.com**.

