Infrared Therapy Sauna Blanket

**Multilingual
FR, IT,
DE & ES**

Product User Manual

www.nebula-appliances.com

Elevate Your Wellness Experience

Welcome to the world of ultimate relaxation with the Infrared Sauna Blanket! Designed with advanced technology and ancient healing principles, this product elevates your wellness experience. In this manual, find setup, usage, and maintenance instructions. Immerse yourself in the soothing warmth and therapeutic benefits of infrared therapy combined with the comforting embrace of a sauna blanket.

IMPORTANT - READ BEFORE USE!

Consult with a medical practitioner before use if you have a medical/other condition listed below, as the use of the sauna blanket is not recommended without your doctor's approval.

- Skin diseases, burns, or reddening
- Sensitivity to heat/temperature
- Heart disease or illness
- Broken bones
- Just had surgery
- Feel nauseous, have a headache or other symptoms
- Splanchnic diseases
- Osteoporosis
- Suffer from abnormal blood pressure
- Currently pregnant or breastfeeding
- Anhidrosis
- Cancer or other malignant conditions
- You have any metal or other body implants

Travel bag included

! Electrical Safety

- Basic safety precautions should be followed to reduce the risk of fire, electric shock, and injury. Misuse or failure to follow directions may cause overheating, fire, or injury to user and/or property.
- Before using the blanket, make sure that your local electric network complies with the voltage and frequency indicated in this manual.
- Otherwise, the machine may be damaged in case of excessive current or voltage rise.
- It is dangerous to change or try to modify the function of this product in any way.
- **Do NOT** attempt to operate this unit if the power cord is frayed or broken.
- **Do NOT** spill any liquids into or onto your unit.
- During long periods of non-use, disconnect the unit's main power.
- Controller is not waterproof, **Do NOT** immerse in any liquids and always keep away from wet or damp areas.
- **Do NOT** use close to other electrical heated products or in high temperature places.

Safety Instructions

- **Do NOT** wear metal jewellery or any other metal items when using the blanket.
- **Do NOT** put the blanket on a leather bed, please put on a cotton towel between the blanket and bed. The high temperatures can cause the PVC and leather to fuse together.
- **Do NOT** use corrosive cleaning products to clean the blanket.
- **Do NOT** let the blanket heat without a person inside the blanket. This can damage the blanket.
- **Do NOT** fold the blanket when it is still hot. Allow the blanket to cool completely before storing.
- **Do NOT** wrap the cord tightly around the mat or controller. Allow the cord to be loose.
- **Do NOT** store items on the mat, this could potentially damage the blanket.
- **Do NOT** put the blanket in direct sunlight. When storing, store the blanket in a cool, dark place.
- Please keep the inside of the blanket clean and dry.
- Make sure you can easily unfasten the velcro on the blanket before turning the power on.
- Remove the power cord from the control box when not in use. **Do NOT** leave the working sauna unattended.
- **Do NOT** allow children, toddlers, or infants to use this blanket.

First Time Use

What do I wear?

- It is recommended to wear comfortable, casual, skin-covering clothing.
- We suggest breathable attire.
- Your feet also need to be covered, so socks must be worn.


Where do I use it?

- The sauna blanket can be placed on any heat-resistant surface like your bed or sofa (No waterbeds, latex, plastic, memory foam, mechanically adjusted beds, bunk beds or child's bed/crib).
- **Tip** – We recommend putting a yoga/fitness mat underneath the blanket, and make sure you are comfortable by adding a pillow/cushion. A suggestion could be the **Nebula Contoured Orthopaedic Pillow (Available on Amazon)**.
- Before each use, inspect the mat for any damage. Inspect the controller and wires to ensure there is no fraying and do **NOT** use if found.

CAUTION: Do NOT exceed 50 minutes in the sauna blanket. **Do NOT** use this heated product directly on synthetic, latex, or flammable items. If you need to stop the session at any time, just press the 'Power' key.

1. Place the sauna blanket on a bed or other flat area. **Do NOT** use the blanket in a wrinkled condition.
2. Connect the sauna blanket to the controller and then connect the controller to the power supply. Check to ensure both connections are secure.
3. Press the power button to activate the controller. **CAUTION: Controller may be hot, do NOT place on top of the blanket.**
4. Select your temperature by hitting the temperature icon button. **Important:** To begin, start at a lower temperature like 60-70 degrees Celsius. The temperature can be increased to a maximum of 80 degrees Celsius. As you become more used to the blanket, then raise the temperature to find the optimal level for you.
5. Set the desired amount of time you want to use the blanket for. Do not exceed more than 50 minutes of use in the blanket.
 - I. The controller will automatically shut off the blanket after 60 minutes of use.
 - II. If at any time during use, you feel overheated or uncomfortable before the time is up, please exit the mat.
 - III. Disconnect from the power supply.
6. Press the start button to begin the temperature and timing settings.

Important: If you begin using high temperatures and feel too hot, lower the temperature. If you need to stop the session early, simply exit the blanket and switch the power button off.

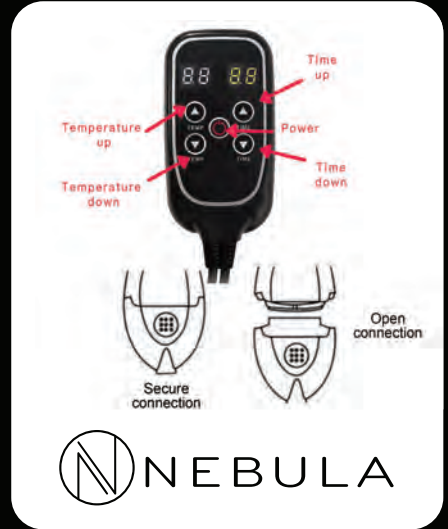Cleaning and Storage Instructions

- The blanket is waterproof, but we do suggest cleaning in between each use to keep the blanket in great condition.
- Use a soft cloth or towel to absorb any sweat left behind.
- Use a non-toxic cleaner (gentle soap) or non-toxic cleaning wipe.
- **Do NOT** use a strong corrosive cleaner or alcohol to clean the sauna blanket.

Suggested Storage Tips

- Please ensure the sauna blanket has completely cooled down before cleaning and storing it.
- Store in a cool, dry place.
- **Do NOT** place any objects on top of the sauna blanket.
- **Do NOT** submerge the sauna blanket in water.
- **Do NOT** place in a washing machine or dry clean.

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